



BIOCHEMICAL OXYGEN DEMAND (BOD) SOLUTIONS

BOD PRO SYSTEM FEATURES

The MANTECH BOD Pro™ system automates the 5-day and 7-day Biological Oxygen Demand (BOD) analysis methods approved for regulatory compliance around the world.

BOD PRO™ CONFORMS TO THE FOLLOWING REGULATORY NORMS:

- ASTM Standard Methods 5210B, 23rd edition
- EPA 405.1
- ISO EN 5815-1
- ISO EN 1899-1, -2



- 1 AM400™ Autosampler
 Maintenance- and binding-free autosamplers with varying capacities can fit all system components on platform.
- 2 YSI MultiLab DO Meter
 Single, dual or triple membrane or optical DO probes.
- 3 Sample Preparation
 Precise dilution water, seed and/or inhibitor (ATU) addition.
- 4 IntelliRinse™


 Dynamic rinse station to clean probe(s) in between each dilution or sample set.
- Racks
 12 or 18 position racks to accommodate any bottle size.
- 6 BOD Pro™ Software

 MANTECH coded software with intuitive user interface.

AUTOMAX 400 SERIES

MANTECH's Next Generation Autosamplers



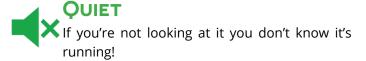






with a single anchor point to eliminate binding.

ONE SIZE DOES NOT FIT ALL

Autosamplers range from 18 - 108 300mL bottle capacities with unlimited batch sizes - BOD Pro software will prompt analysts via email notification to swap out racks - with 12 or 18 position racks options available.

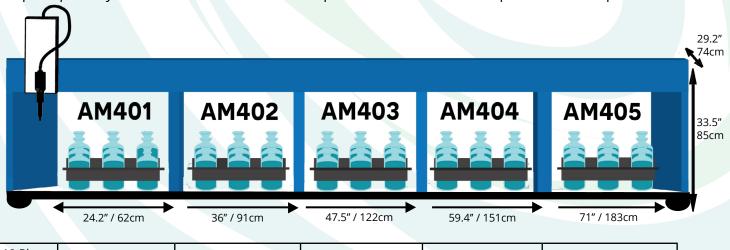

12 Place	N/A	36	60	84	108
Racks		(3 Racks)	(5 Racks)	(7 Racks)	(9 Racks)
18 Place	18	36	54	72	90
Racks	(1 Rack)	(2 Racks)	(3 Racks)	(4 Racks)	(5 Racks)

Customizable racks options available to accommodate any bottle size.

"...built like a tank!"

- A MANTECH Customer

SELECT YOUR AUTOMATION LEVEL

MANTECH systems can be configured to automate as much or as little of the biochemical oxygen demand analysis freeing up analyst's time and generating efficiencies for your laboratory.

SAMPLE PREPARATION

MANTECH's peristaltic pump, IntelliPump™, can automate dilution water, seed and/or inhibitor (ATU) addition.

- Software-controlled calibrations and flow rate/direction
- Precision +/-1%
- Plug-and-play functionality
- Flow rate: 12.5 to 860 mL/min (clockwise & counterclockwise)

Looking for higher precision?

100,000 step resolution syringebased pump can deliver precision +/-0.01%!





PH ADJUSTMENT

As per SM 5210B 23rd edition, if pH is not between 6.0 - 8.0, then sample must be adjusted to pH between 6.5 - 7.5. This system measures and adjusts the pH of the sample while recording all data.

pH METER



pH Meter w/ Adjustment



AUTOMATED pH & ADJUSTMENT



SAMPLE TRANSFER

MANTECH automates pipetting field samples into dilution bottles with IntelliDose™, a patented piston-based pump to deliver precise aliquots without clogging.

- Add automated pH adjustment
- No syringes, valves, sample loops or levelling required
- Flexible configuration
- Use autosampler for additional analysis (pH, alk, EC, etc.) during downtime
- Customizable racks for field bottles





BOD ANALYSIS

BOD PRO™ SOFTWARE

The MANTECH designed and coded BOD Pro™ software is easy to use and is designed to comply with user protocol requirements, while allowing user-defined flexibility.



ENSURES COMPLIANCE

- Tracks incubations periods
- Performs all calculations with ability to modify equations and methods

AGNOSTIC TO ANY LIMS

- Full import and export capabilities
- Report editor add 10 user-defined data columns

MULTILINGUAL FUNCTIONALITY

Alternate between any language

USER-FRIENDLY INTERFACE

- Shortcut buttons for routine runs, PM and more
- Visual rack layout
- Use familiar Excel functionality for run set-up

SAFEGUARD YOUR DATA

- Designed for WIN 10 & 11 Pro OS
- Account permissions for each analyst

MULTITASKING CAPABILITY

Access all functionality during a run



Manually prepare dilutions and place DO probe in bottles and BOD Pro™ software will take care of the rest!

- Reporting
- Calculations
- Bottle/incubation tracking

HIGHLY CONFIGURABLE

MANTECH systems are designed to grow with your laboratory's changing demands with additional features and automation to be upgraded at any time.

MULTIPLE PROBES

Single, dual or triple DO probe options:

- Up to x3 bottles at a time with sample preparation
- Up to x6 bottles at a time with DO readonly system



QC RACK

Keep your racks available for more dilutions with a separate QC rack for:

- GGA's
- Blanks
- Probe park



BARCODE SCANNER

Integrate a barcode scanner to BOD Pro software to keep track of your dilutions.



10-MINUTE BOD ESTIMATOR

MANTECH's patented BOD estimator delivers results in 10-minutes or less which enables operators to:

- Reduce the number of dilutions for samples with little to no historical data
- Get results back to engineers to make informed treatment decisions with confidence
- Ensure valid final BOD5 results for all samples



WHY MANTECH

TRUSTED WORLDWIDE



With 3,000+ analyzers in the market, MAN-TECH systems are trusted to produce more than 670,000 results every day that protect public health and the environment.

SUPPORT

MANTECH offers several options:

- Context sensitive GIFs and instructions built into software
- On-Demand virtual support
- In-person visits with MANTECH Technical Representatives or distributors

END-TO-END MANUFACTURER



MANTECH develops and manufactures all hardware and software as a single source provider.

IOT CONNECTIVITY



All MANTECH system components feature IoT connectivity via USB or ethernet/IP connection for bi-directional communication with the software.

DUAL TITRATION CAPABILITY

MANTECH's BOD systems can be configured to add titration and multi-parameter functionality (pH, EC, alkalinity, etc.) in the same platform.

LEARN MORE

MODULAR DESIGN



The modular design allows for MANTECH systems to be highly configurable and plugand-play for system upgrades and trouble-shooting.

5473 HIGHWAY 6 NORTH GUELPH, ONTARIO, CANADA N1H 6J2

WWW.MANTECH-INC.COM INFO@MANTECH-INC.COM

- +1 (519) 763-4245
- +1 (866) 763-2122

@mantech inc