



Solutions for Environmental Testing

Soil Application — In Accordance with US EPA 3545A or US EPA 3546



EDGE and MARS 6

Pressurized Fluid Extraction

Microwave Solvent Extraction

We Simplify Science



Options to fit your needs.

EDGE[®]



MARS 6[™]



Sequential	Batch
Mixed Sample Types	Like Sample Types
No Solvent Limitations	Solvent Must Have Dipole Moment
Filtered Samples	Samples Must be Manually Filtered
Extraction Only	Extraction and Digestion in One System
Great For One-Off Difficult Samples	Great For High Throughput




EDGE

Automated Solvent Extraction System

The EDGE is a sequential automated pressurized fluid extraction system that is faster than Soxhlet and simpler than other automated extraction systems. With the Q-Cup® Sample Holder, it is perfect for efficiently preparing samples for US EPA 3545A. The EDGE automates the extraction and washing process and provides a filtered and cooled extract.

Sample Holder
Q-Cup

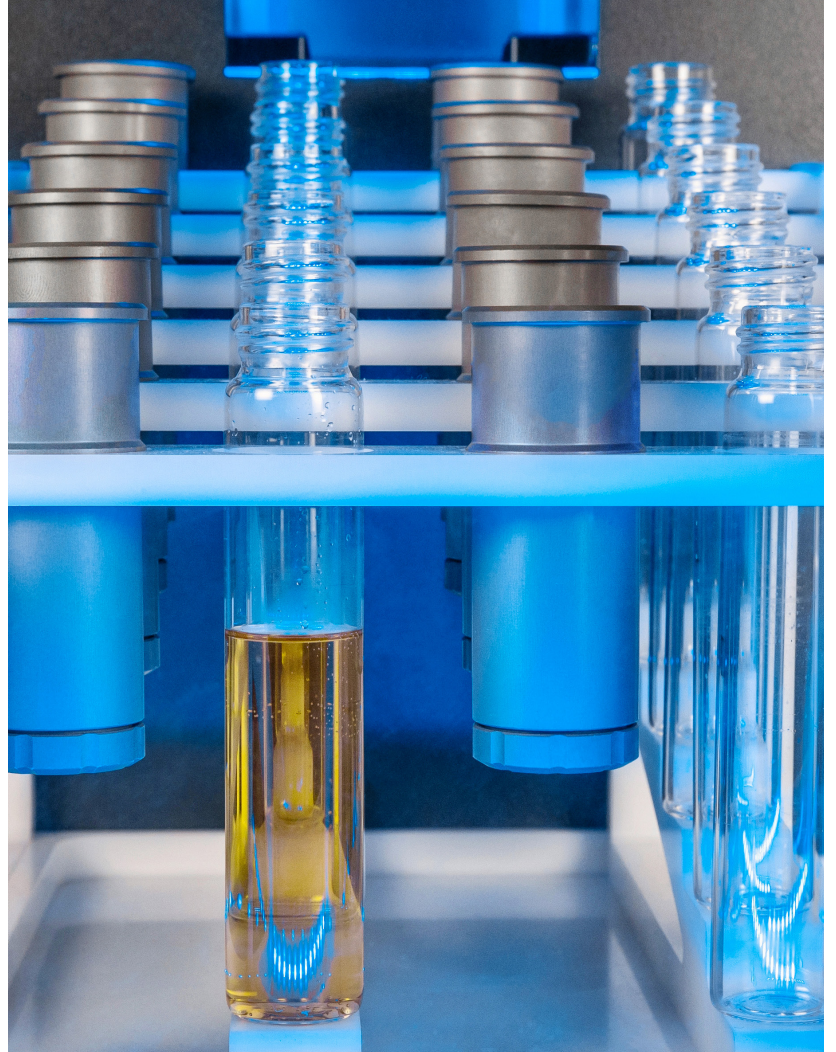



MARS 6

Microwave Extraction System

The MARS 6 is a batch microwave extraction system with an intelligent design that's remarkably easy to use. With an industry leading 2000 watts of installed power, it is the perfect tool for preparing samples for US EPA 3546. The MARS 6 can easily prepare up to 40 samples in a single batch, or up to 24 samples using disposable glass inserts, providing the highest throughput of any microwave extraction system on the market.

Vessels


MARSXpress Plus



Results

Both the EDGE and the MARS 6 extracted 20 g of a BNA CRM. Acceptable recoveries and RSD values were obtained.

MARS 6 and EDGE CRM Recovery

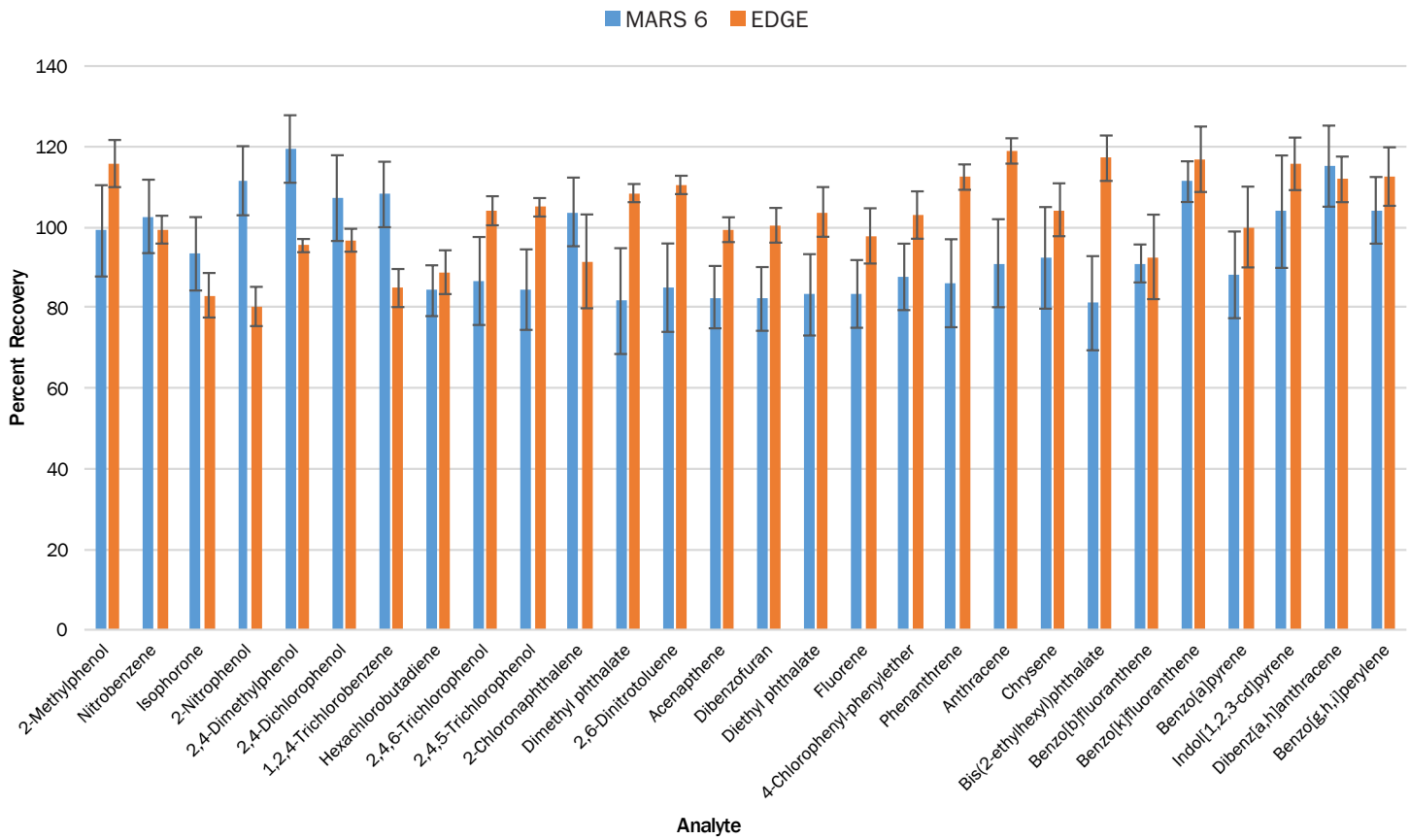





We Simplify Science

cem.com



United States Headquarters: +1 (704) 821-7015 | info@cem.com
For distributors and subsidiaries in other regions, visit cem.com/contact